

# PHILIPS

## HearLink

In-the-Ear  
Hearing Aids



Made for

 iPhone | iPad | iPod

Works with

android 

# Instructions for Use

## Model overview

This booklet is valid for the HearLink family in the following hearing aid models, battery sizes, and styles:

### FW 1.4

- HearLink 9030      GTIN: (01) 05714464059997  
(HEI9030, HEI9031, HEI9032, HEI9033, HEI9034)
- HearLink 7030      GTIN: (01) 05714464060009  
(HEI7030, HEI7031, HEI7032, HEI7033, HEI7034)
- HearLink 5030      GTIN: (01) 05714464060016  
(HEI5030, HEI5031, HEI5032, HEI5033, HEI5034)
- HearLink 3030      GTIN: (01) 05714464060023  
(HEI3030, HEI3031, HEI3032, HEI3033, HEI3034)
- HearLink 2030      GTIN: (01) 05714464060030  
(HEI2030, HEI2031, HEI2032, HEI2033, HEI2034)

## Battery size

- 10    312

**Bluetooth®**    Yes    No

## Styles

- |                                 |                         |
|---------------------------------|-------------------------|
| <input type="checkbox"/> IIC    | Invisible-In-the-Canal  |
| <input type="checkbox"/> CIC    | Completely-In-the-Canal |
| <input type="checkbox"/> ITC    | In-the-Canal            |
| <input type="checkbox"/> ITE HS | In-the-Ear Half Shell   |
| <input type="checkbox"/> ITE FS | In-the-Ear Full Shell   |

## Introduction to this booklet

This booklet guides you on how to use and maintain your new hearing aids. Ensure you read the booklet carefully including the **Warnings** section. This will help you get the most benefit of your new hearing aid.

More information is available on [www.hearingsolutions.philips.com](http://www.hearingsolutions.philips.com)

Your hearing care professional has adjusted the hearing aid to meet your needs. If you have additional questions, contact your hearing care professional.

A hearing care professional\* (hearing aid professional, audiologist, ENT (ear, nose and throat) doctor, and hearing aid dispenser) is a person who is appropriately educated and has proven competency in professionally assessing hearing, selecting, fitting, and delivering hearing instruments and rehabilitation care to persons with hearing loss.

The education of the hearing care professional is in accordance with national or regional regulations.

\* The job title may vary from country to country.

For your convenience this booklet contains a navigation bar to help you navigate easily through the different sections.

| **About** | [Getting started](#) | [Daily use](#) | [Options](#) | [Warnings](#) | [More info](#) |

## Intended use

<b>Intended use</b>	The hearing aid is intended to amplify and transmit sound to the ear.
<b>Indications for use</b>	Bilateral or unilateral impaired hearing of sensorineural, conductive, or mixed type ranging from a slight (25 dB HL*) to severe (80 dB HL*) degree of hearing loss, with an individual frequency configuration.
<b>Intended user</b>	Person with hearing loss using a hearing aid and their caregivers. Hearing care professional responsible for adjusting the hearing aid.
<b>Intended user group</b>	Adults.
<b>Use environment</b>	Indoor and outdoor.
<b>Contraindications</b>	Not suitable for infants below 36 months.** Users of active implants must pay special attention when using the hearing aid. For more information read the <b>Warnings</b> section.
<b>Clinical benefits</b>	The hearing aid is designed to provide better speech understanding to help ease communication with the aim of improving quality of life.

\* As specified by the American Speech-Language-Hearing Association, *asha.org*, using pure-tone average of 0.5, 1 and 2 kHz.

\*\* In accordance with current Pediatric Amplification guidelines, custom hearing aids are not typically considered for children of any age. However, a hearing care professional may, in special circumstance and on a case-by-case basis, select them for children older than 36 months.

**IMPORTANT NOTICE**

The hearing aid amplification is uniquely adjusted and optimised to your personal hearing capabilities during the hearing aid fitting performed by your hearing care professional.

# Table of contents

## About

Identify your hearing aid	12
Hearing aids with size 10 battery	13
Hearing aids with size 312 battery	14
Identify left and right hearing aid	15
MultiTool for handling batteries and cleaning	16

## Getting started

Turn the hearing aid ON and OFF	18
Replace the battery	19
How to replace the battery	20

## Daily use

Insert the hearing aid	22
Remove your hearing aid	23
Caring for your hearing aid	24
Filter replacement	26
Replace ProWax miniFit filter	27



Replace T-Cap filter	28
Replace O-Cap filter	29
Change programs	30
Change volume with push-button	32
Change volume with volume wheel	33

## **Options**

Mute the hearing aid	34
Flight mode	35
Use hearing aids with iPhone, iPad, and iPod touch	36
Use your hearing aids with Android devices	37
Wireless accessories	38
Other options	40
Optional features and accessories	41

## **Warnings**

General warnings	42
------------------	----

# Table of contents

## More info

Troubleshooting guide	48
Water & dust resistant (IP68)	50
Conditions of use	51
Warranty	52
International warranty	53
Technical information	54
Your individual hearing aid settings	59




## Identify your hearing aid

For your in-the-ear hearing aid, there are several different styles available with different battery sizes. Please see the model overview for your style and battery size. This will make it easier for you to navigate through this booklet.

# Hearing aids with size 10 battery

(CIC shown)

## What it is and does

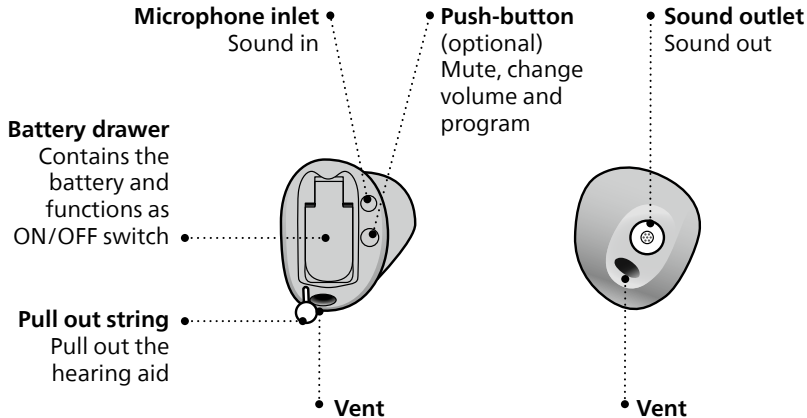


Components may be positioned differently on your hearing aid.

# Hearing aids with size 312 battery


(ITC shown)

## What it is and does

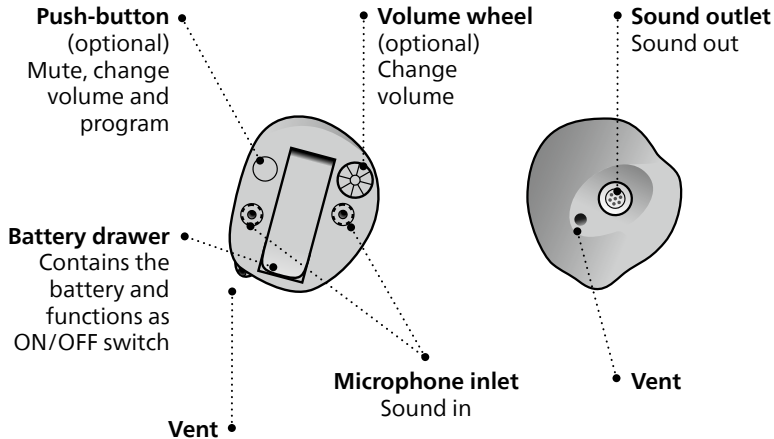


## Identify left and right hearing aid

It is important to distinguish between the left and the right hearing aid as they might be shaped and programmed differently.

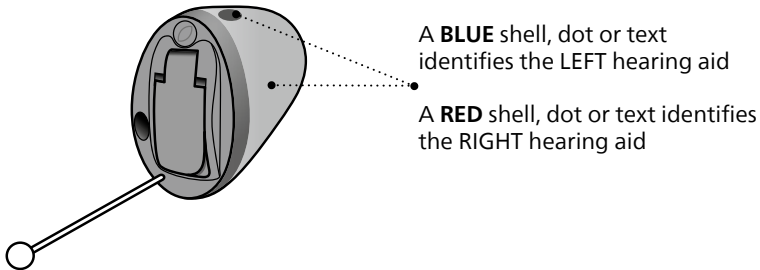

## MultiTool for handling batteries and cleaning

The MultiTool contains a magnet that makes it easier to replace the battery in the hearing aid. It also contains a brush and a wire loop for cleaning and removing earwax. If the vent is very small, a special tool may be required. If you need a new tool please contact your hearing care professional.

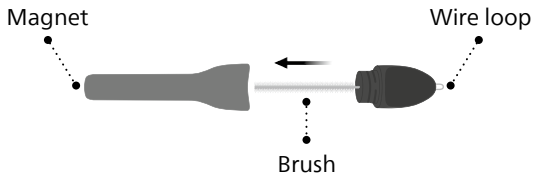

### IMPORTANT NOTICE

The MultiTool has a built-in magnet. Keep the MultiTool at least 30 centimetres (1 foot) away from credit cards and other magnetically sensitive devices.






## Turn the hearing aid ON and OFF

The battery drawer is also used to switch the hearing aid ON and OFF. To save battery life, make sure your hearing aid is switched OFF when you are not wearing it. If you wish to return to the standard settings of the hearing aid, simply open and then close the battery drawer (quick reset).


### Turn ON

Close the battery drawer with the battery in place.

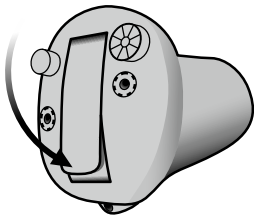

### Turn OFF

Open the battery drawer.

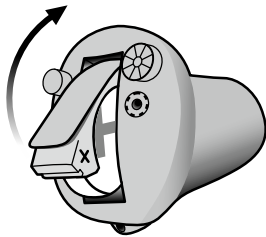


**Note:** Open the battery door fully to make sure the hearing aid does not drain the battery fully and to allow air to circulate whenever you are not using your hearing instrument, especially at night or for longer periods of time.

## Replace the battery

When it is time to replace the battery, you will hear three beeps repeated at moderate intervals until the battery runs out.

 **Three beeps\***  
= The battery is running low

 **Four beeps**  
= The battery has run out

### Battery tip


To make sure the hearing aid is always working, bring spare batteries with you, or replace the battery before you leave home.

\* If your hearing aid has Bluetooth® this will be turned off and it will not be possible to use wireless accessories.

**Note:** Batteries need to be replaced more often if you are streaming audio or music to your hearing aids.


# How to replace the battery

## 1. Remove

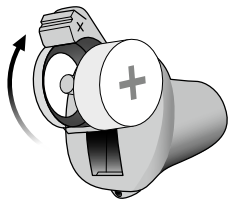

Fully open the battery drawer. Remove the battery.


## 2. Uncover

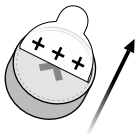

Remove the sticky label from the + side of the new battery.


Tip:  
Wait 2 minutes so that the battery can draw air, to ensure optimal functioning.

## 3. Insert

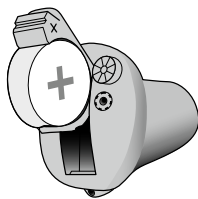


Insert the new battery into the battery drawer. Ensure that the + side of the battery and battery drawer align.

#### 4. Close

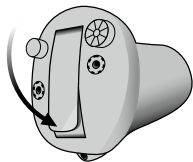

Close the battery drawer. The hearing aid will play a jingle. Hold the hearing aid close to your ear to hear the jingle.

#### Tip

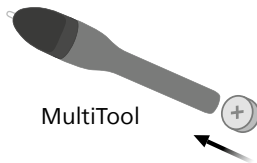

The MultiTool can be used for battery change. Use the magnetic end to remove and insert batteries.


# Insert the hearing aid

## Step 1



Place the tip of the hearing aid in your ear canal.

## Step 2

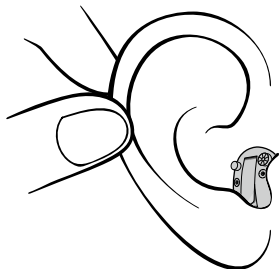

Gently pull your ear outwards and push the hearing aid into the ear canal, twisting slightly if necessary. Follow the natural contour of the ear canal. Push the hearing aid to make sure it fits comfortably in the ear.

## Remove your hearing aid

Hold the hearing aid by the pull-out string (if available). Gently pull the hearing aid from the ear canal.

If your hearing aid doesn't have a pull-out string, you can remove it by pulling on the edge of the hearing aid.

### **IMPORTANT NOTICE**


DO NOT use the battery door as a handle to insert or remove your hearing aids, as it is not designed for this purpose.

## Caring for your hearing aid

When handling your hearing aid, hold it over a soft surface to avoid damage if you drop it.


### Cleaning the hearing aid

Carefully brush away debris from the microphone inlets with a clean brush. Gently brush the surface. Make sure that the filters do not fall off.

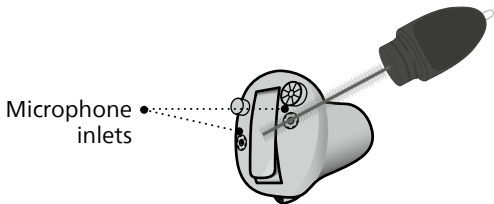




Clean the vent by pressing the brush through the hole while twisting it slightly.

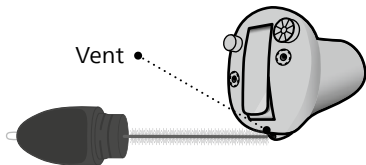

- If the vent is very small, a special tool may be required. Please consult your hearing care professional.

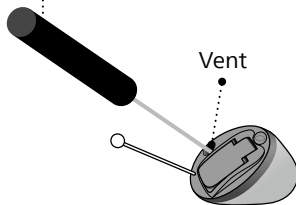

### **IMPORTANT NOTICE**

To clean the hearing aids, use a soft, dry cloth. The hearing aids must never be washed or immersed in water or other liquids.

## Filter replacement

The filters keep wax and debris from damaging the hearing aid. If the filters become clogged, please replace the filters or contact your hearing care professional.

- ProWax miniFit filter protects the sound outlet.
- O-cap and T-cap protect the microphone inlet.

Please refer to the following pages for instructions on how to replace the appropriate filters.


### **IMPORTANT NOTICE**

Ensure you always use the same type of wax filter as was originally supplied with the hearing aids.

If you are in doubt about the use or replacement of ProWax miniFit filters, contact your hearing care professional.


# Replace ProWax miniFit filter

## 1. Tool

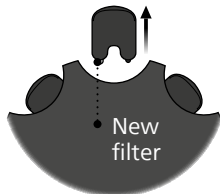


Remove the tool from the shell. The tool has two pins, one empty for removal and one with the new ProWax miniFit filter.

## 2. Remove

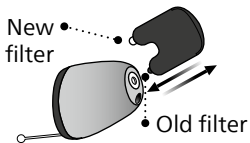

Push the empty pin into the ProWax miniFit filter in the hearing aid and pull it out.

## 3. Insert

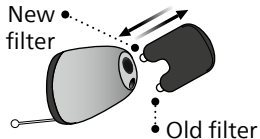


Insert the new ProWax miniFit filter using the other pin, remove the tool and throw it out.

# Replace T-Cap filter


(hearing aids with size 10 battery)

## 1. Tool

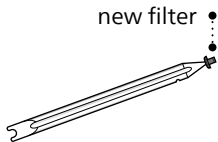


Remove the tool from the packaging. The tool has two ends, one for removal and one with the new T-Cap filter.

## 2. Remove

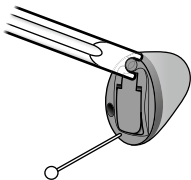

Push the tool fork under the top edge of the used T-Cap filter and lift it out.

## 3. Insert

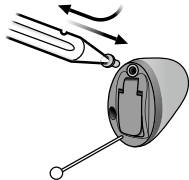


Insert the new T-Cap filter and remove the tool by twisting it slightly. Throw the tool out after use.

# Replace O-Cap filter


(hearing aids with size 312 batteries)

## 1. Tool

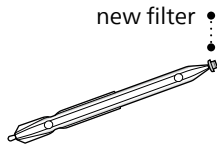


Remove the tool from the packaging. The tool has two ends, one for removal and one with the new O-Cap filter.

## 2. Remove

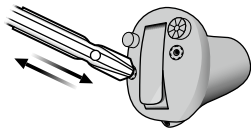

Push the pointed end of the tool into the existing O-Cap filter and pull it out.


## 3. Insert

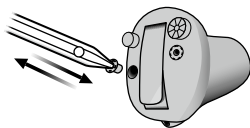

Insert the new O-Cap filter using the other end of the tool and remove. Throw the tool out after use.

## Change programs

Your hearing aid can have up to 4 different programs. These are programmed by your hearing care professional.

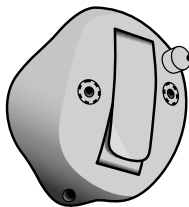





- Press the button to change program. Use a short press if the push-button is used for program change only and a long press if it is also used for volume control.

Note that if you have two synchronised hearing aids (both hearing aids respond when either push-button is operated), the RIGHT hearing aid switches forward from e.g., program 1 to 2 and the LEFT hearing aid switches backwards from e.g., program 4 to 3.

If your hearing aids work independently, you must press the buttons on each hearing aid.


## To be filled out by the hearing care professional

Program	Sound you will hear when activated	When to use
1	 "1 beep"	
2	 "2 beeps"	
3	 "3 beeps"	
4	 "4 beeps"	

Program change:	
<input type="checkbox"/> Independent	<input type="checkbox"/> Synchronised
<input type="checkbox"/> LEFT	<input type="checkbox"/> RIGHT
<input type="checkbox"/> Short press	<input type="checkbox"/> Long press

## Change volume with push-button

The push-button allows you to adjust the volume.  
You may hear a click when you turn the volume up or down.

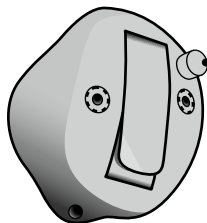

A short press on the  
RIGHT hearing aid  
increases the volume

A short press on the  
LEFT hearing aid  
decreases the volume



### To be filled out by the hearing care professional


Volume change	<input type="checkbox"/> LEFT	<input type="checkbox"/> RIGHT
---------------	-------------------------------	--------------------------------



## Change volume with volume wheel


The volume wheel allows you to adjust the volume.

You may hear a click when you turn the volume up or down.

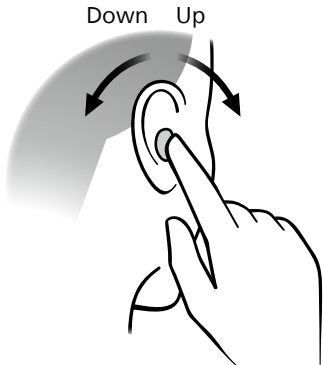


Turn forward to  
turn up volume

Turn backwards  
to turn down the  
volume

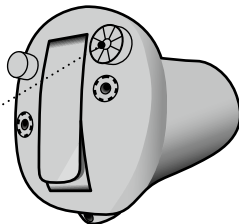

## Mute the hearing aid

Use the mute function if you need to silence the hearing aid. Only available for hearing aids with push-button.

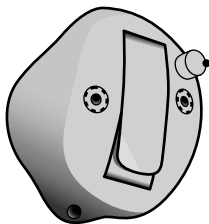

- Apply a very long press (> 2.4 seconds) to the push-button to mute the hearing aid. To reactivate the hearing aid, push the button briefly.

NOTE: the mute function only mutes the microphone(s) on the hearing aids.

### IMPORTANT NOTICE

Do not use the mute function as an off switch, as the hearing aids are still using battery power in this mode.

## Flight mode


(hearing aids with Bluetooth and push-button)

When boarding an airplane or entering an area in which it is prohibited to radiate radio signals, e.g. during flight, flight mode must be activated. The hearing aid will still be working. It is only necessary to activate flight mode on one hearing aid to turn off Bluetooth on both hearing aids. If your hearing aid does not have a push-button you need to turn off your hearing aid.

**To activate and deactivate**

Press the push-button for at least 7 seconds. A jingle confirms your action.

Opening and closing the battery drawer will also deactivate flight mode.

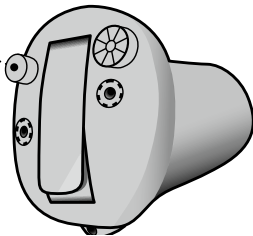

## Use hearing aids with iPhone, iPad, and iPod touch

Your hearing aids are Made for iPhone® and allow for direct communication and control with iPhone, iPad® or iPod touch®.

For assistance in using these products with your hearing aids, contact your hearing care professional or visit:  
[hearingsolutions.philips.com/hearing-aids/connectivity](https://hearingsolutions.philips.com/hearing-aids/connectivity)

To see how to pair your hearing aids with your iPhone visit:  
[hearingsolutions.philips.com/support/how-to-faq/](https://hearingsolutions.philips.com/support/how-to-faq/)

For information on compatibility, visit:  
[hearingsolutions.philips.com/compatibility](https://hearingsolutions.philips.com/compatibility)

Use of the Made for Apple badge means that an accessory has been designed to connect specifically to the Apple products identified in the badge, and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards.

Note that use of this accessory with iPod, iPhone, or iPad may affect wireless performance.

## **Use your hearing aids with Android devices**

Your hearing aids support Audio Streaming for Hearing Aids (ASHA) and allow for direct streaming from selected Android™ devices.

For assistance in using your hearing aids with Android devices, contact your hearing care professional or visit:  
[hearingsolutions.philips.com/hearing-aids/connectivity](https://hearingsolutions.philips.com/hearing-aids/connectivity)

To see how to pair your hearing aids with your Android device visit:  
[hearingsolutions.philips.com/support/how-to-faq/](https://hearingsolutions.philips.com/support/how-to-faq/)

For information on compatibility, visit:  
[hearingsolutions.philips.com/compatibility](https://hearingsolutions.philips.com/compatibility)

## Wireless accessories

As an enhancement to your hearing aids, a range of wireless accessories are available. These can enable you to hear and communicate better in many everyday situations.

### AudioClip

When paired with your mobile phone, you can use the hearing aids as a hands-free headset, or as a remote microphone.

### TV Adapter

Wireless transmitter of sound from TV and electronic audio devices. TV Adapter streams sound directly to your hearing aids.

### Remote Control

Offers the ability to change program, adjust volume, or mute your hearing aids.

### HearLink app

An application that lets you control your hearing aid from your mobile phone or tablet. For iPhone, iPad, iPod touch, and Android™ devices. Ensure that you only download and install the app from the official app stores.

For more information visit:  
[hearingsolutions.philips.com/hearing-aids/connectivity](https://hearingsolutions.philips.com/hearing-aids/connectivity)

Apple, the Apple logo, iPhone, iPad, and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc. Android, Google Play, and the Google Play logo are trademarks of Google LLC.

## Other options

### Telecoil - optional for select styles

Telecoil helps you hear better when using a telephone with a built-in loop or when you are in buildings with teleloop systems such as theatres, churches, or lecture rooms. This symbol or a similar sign is shown wherever a teleloop has been installed.



For more information, please contact your hearing care professional.



## Optional features and accessories

The features and accessories described on the following pages are optional. Please contact your hearing care professional to find out how your hearing aid is programmed.

If you experience difficult listening situations, a special program may be helpful. These are programmed by your hearing care professional.

*Write down any hearing situations in which you may need help.*

---

---

---

## **General warnings**

For your personal safety and to ensure correct usage, you should familiarise yourself fully with the following general warnings before using your hearing aids. Consult your hearing care professional if you experience unexpected operations or serious incidents with your hearing aid during use or because of its use. Your hearing care professional will support you with issue handling and, if relevant, reporting to the manufacturer and/or the national authorities.

Note that hearing aids do not restore normal hearing and do not prevent or improve a hearing impairment resulting from organic conditions.

Hearing aids are only a part of hearing habilitation and may need to be supplemented by auditory training and instruction in lipreading.

Furthermore, note that in most cases, infrequent use of a hearing aids does not permit a user to attain full benefit from it.

### **Usage of hearing aids**

Hearing aids should only be used as directed and adjusted by your hearing care professional. Misuse can result in sudden and permanent hearing loss.

Never allow others to wear your hearing aid, as incorrect usage could cause permanent damage to their hearing.

### **Choking hazards and risk of swallowing batteries or other small parts**

Hearing aids, their parts and batteries should be kept out of reach of children and anyone who might swallow these items or otherwise cause injury to themselves.

Batteries have occasionally been mistaken for pills. Therefore, check your medicine carefully before swallowing any pills.

*If a battery, hearing aid or small part is swallowed, see a doctor immediately.*

Continues on next page

### **Battery use**

Always use batteries recommended by your hearing care professional. Batteries of low quality may leak and cause bodily harm.

Never attempt to recharge your batteries and never dispose of batteries by burning them. There is a risk that the batteries will explode.

### **Dysfunction**

Be aware of the possibility that your hearing aids may stop working without notice. Keep this in mind when you depend on warning sounds (e.g. when you are in traffic). The hearing aids may stop functioning, for instance if the batteries have expired or if the tubing is blocked by moisture or earwax.

### **Active implants**

The hearing aid has been thoroughly tested and characterised for human health according to international standards for human exposure (Specific Absorption Ratio - SAR), induced electromagnetic power and voltages into the human body.

The exposure values are well below international accepted safety limits for SAR, induced electromagnetic power and voltages into the human body defined in the standards for human health and coexistence with active medical implants, such as pacemakers and heart defibrillators.

The Autophone magnet or MultiTool (which has a built-in magnet) should be kept more than 30 centimetres (1 foot) away from the implant, e.g. do not carry it in your breast pocket. Follow the guidelines recommended by the manufacturers of implantable defibrillators and pacemakers on use with magnets.

## **General warnings**

If you have an active brain implant, contact the manufacturer of your implantable device for information about the risk of disturbance.

### **X-ray/CT/MR/PET scanning, electrotherapy and surgery**

Remove your hearing aid before X-ray examinations and CT/MR/PET scans, electrotherapy, surgery, etc. as your hearing aid may be damaged when exposed to strong electromagnetic fields.

### **Heat and chemicals**

The hearing aid must never be exposed to extreme heat, e.g. left inside a parked car in the sun.

The hearing aid must not be dried in microwave ovens or other ovens.

The chemicals in cosmetics, hairspray, perfume, aftershave lotion, suntan lotion, and insect repellent can damage the hearing aid. Always remove your hearing aid before applying such products and allow time to dry before use.

### **Power hearing aid**

Special care should be exercised in selecting, fitting, and using a hearing aids where the maximum sound pressure capability exceeds 132 dB SPL (IEC 60318-4 / IEC 711), as there may be a risk of impairing the remaining hearing of the hearing aid user. For information on whether your instrument is a power hearing aid, see the **Your individual hearing aids settings** section in this booklet.

Continues on next page

### **Possible side effects**

Hearing aids and earpieces may cause an accelerated accumulation of earwax.

The non-allergenic materials used in hearing aids may, in rare cases cause a skin irritation or other side effects.

*If these conditions occur, seek consultation with a physician .*

### **Use on aircraft**

If your hearing aids have Bluetooth wireless technology, when on board an aircraft, the hearing aids must be put into flight mode to deactivate Bluetooth, unless Bluetooth is permitted by the flight personnel.

The hearing aid needs to have a push-button for you to be able to activate flight mode, otherwise you must turn off your hearing aids.

*To find out whether your hearing aids have Bluetooth wireless technology, see the model overview in your hearing aids' Instructions for use.*

### **Connection to external equipment**

The safety of the hearing aids, when connected to external equipment (with an auxiliary input cable and/or with USB cable and/or directly), is determined by the external signal source. When the hearing aid is connected to external equipment which is plugged into a power socket, this equipment must comply with IEC 62368-1 or equivalent safety standards.

### **Use of third-party accessories**

Only use accessories, transducers, or cables supplied by the manufacturer. Non-original accessories may result in reduced electromagnetic compatibility (EMC) of your hearing aid.

## **General warnings**

### **Modification of hearing aids is not allowed**

Changes or modifications not expressly approved by the manufacturer will void the warranty of the equipment.

### **Interference**

The hearing aids have been thoroughly tested for interference according to the most stringent international standards.

Electromagnetic interference may occur in the vicinity of equipment with the symbol on the left. Portable and mobile RF (radio frequency) communications equipment can affect the performance of the hearing aids. If your hearing aids are affected by electromagnetic interference, move away from the source to reduce the interference.



# Troubleshooting guide

Symptom	Possible causes
<b>No sound</b>	Worn-out battery
	Clogged sound outlet
	Clogged microphone inlet
	Hearing aid microphone muted
<b>Intermittent or reduced sound</b>	Clogged sound outlet
	Moisture
	Worn-out battery
<b>Squealing noise</b>	Hearing aid not inserted properly
	Ear wax accumulated in ear canal
<b>Pairing issue with Apple device</b>	Bluetooth connection failed
	Only one hearing aid paired



## Solutions

Replace the battery

Clean sound outlet\* or replace the wax filter

Clean microphone inlet\* or replace filter (T-Cap or O-Cap)

Unmute the hearing aid microphone

Clean sound outlet or replace the wax filter

Wipe battery and hearing aid with a dry cloth

Replace the battery

Re-insert the hearing aid

Have ear canal examined by your doctor

**1)** Unpair your hearing aid. **2)** Turn Bluetooth off and on again. **3)** Open and close battery drawer on hearing aid. **4)** Re-pair hearing aid (visit [hearingsolutions.philips.com/support/how-to-faq/](https://hearingsolutions.philips.com/support/how-to-faq/))

If none of the above solutions work, consult your hearing care professional for assistance.

\* According to guideline in this booklet

## Water & dust resistant (IP68)

Your hearing aid is dust-tight and protected against ingress of water which means it is designed to be worn in all daily life situations. Therefore you do not have to worry about sweat or getting wet in the rain. Should your hearing aid come in contact with water and stop working, please follow these guidelines:

1. Gently wipe off any water
2. Open the battery drawer and remove the battery and gently wipe off any water in the battery drawer
3. Let the hearing aid dry with the battery drawer left open for approximately 30 minutes
4. Insert a new battery

### **IMPORTANT NOTICE**

Do not wear your hearing aid while showering or participating in water activities. Do not immerse your hearing aid in water or other liquids.

## Conditions of use

<b>Operating Conditions</b>	Temperature: +1°C to +40°C (34°F to 104°F) Humidity: 5% - 93% relative humidity, non-condensing Atmospheric pressure: 700 hPa to 1060 hPa
<b>Transportation and storage conditions</b>	Temperature and humidity shall not exceed the below limits for extended periods during transportation and storage:  Temperature: -25°C to + 60°C (-13°F to 140°F) Humidity: 5% - 93% relative humidity, non-condensing Atmospheric pressure: 700 hPa to 1060 hPa

# Warranty

## Certificate

Name of owner: \_\_\_\_\_

Hearing care professional: \_\_\_\_\_

Hearing care professional's address: \_\_\_\_\_

Hearing care professional's phone: \_\_\_\_\_

Purchase date: \_\_\_\_\_

Warranty period: \_\_\_\_\_ Month: \_\_\_\_\_

Model left: \_\_\_\_\_ Serial no.: \_\_\_\_\_

Model right: \_\_\_\_\_ Serial no.: \_\_\_\_\_

Battery size: \_\_\_\_\_

## International warranty

Your device is covered by an international warranty issued by the manufacturer. This international warranty covers manufacturing and material defects in the device itself, but not in accessories such as batteries, tubing, speakers, earpieces and filters, etc. Problems arising from improper/incorrect handling or care, excessive use, accidents, repairs made by an unauthorised party, exposure to corrosive conditions, physical changes in your ear, damage due to foreign objects entering the device, or incorrect adjustments are NOT covered by the international warranty and may void it. The above international warranty does not affect any legal rights that you might have under applicable national legislation governing the sale of consumer goods in the country where you have bought your device. Your hearing

care professional may also have issued a warranty that goes beyond the clauses of this international warranty. Please consult him/her for further information.

### **If you need service**

Take your device to your hearing care professional, who may be able to sort out minor problems and adjustments immediately. Your hearing care professional may charge a fee for their services.

## Technical information

The hearing aids contain the following two radio technologies:

The hearing aids contain a radio transceiver using short range magnetic induction technology operating at 3.84 MHz. The magnetic field strength of the transmitter is very weak and always below 15 nW (typically below  $-40 \text{ dB}\mu\text{A}/\text{m}$  at a distance of 10 metres ( $-12.20 \text{ dB}\mu\text{A}/\text{ft}$  at a distance of 33 feet).

The hearing aids also contain a radio transceiver using Bluetooth Low Energy and a proprietary, short range radio technology both operating at ISM band 2.4 GHz. The radio transmitter is weak and always below 3 mW equal to 4.8 dBm in total radiated power.

Only use your hearing aids in areas where wireless transmission is permitted.

The hearing aids comply with international standards concerning radio transmitters, electromagnetic compatibility, and human exposure. Due to the limited space available on the hearing aid, relevant approval markings can be found in this booklet. Additional information can be found in the "Specification Guide" on [hearingsolutions.philips.com](http://hearingsolutions.philips.com)

The manufacturer declares that this hearing aid is in compliance with the essential requirements and other relevant provisions of Directive 2014/53/EU.

This medical device complies with Medical Device Regulation (EU) 2017/745.

Declaration of Conformity is available at the headquarters.

SBO Hearing A/S  
Kongebakken 9  
DK-2765 Smørum  
Denmark  
[hearingsolutions.philips.com/doc](https://hearingsolutions.philips.com/doc)

CE 0123



SBO Hearing A/S  
Kongebakken 9  
DK-2765 Smørum  
Denmark



Waste from electronic equipment must be handled according to local regulations.



 **Bluetooth®**

**IP68**

## Description of symbols and abbreviations used in this booklet



### Warnings

Text marked with a warning symbol must be read before using the device.



### Manufacturer

The device is produced by the manufacturer whose name and address are stated next to the symbol. Indicates the medical device manufacturer, as defined in EU Regulation 2017/745.



### CE mark

The device complies with all required EU regulations and directives. The four digit number indicates the identification of the notified body.



### Electronic waste (WEEE)

Recycle hearing aids, accessories or batteries according to local regulations. Hearing aid users can also return electronic waste to their hearing care professional for disposal. Electronic equipment covered by Directive 2012/19/EU on waste and electrical equipment (WEEE).



### Regulatory Compliance Mark (RCM)

The device complies with electrical safety, EMC and radio spectrum requirements for devices supplied to the Australian or New Zealand markets.

IP68

### IP code

Class of protections against harmful ingress of water and particulate matter according to EN 60529. IP6X indicates total dust protection. IPX8 indicates the protection against the effects of continuous immersion in water.



**Bluetooth logo**

Registered trademark of Bluetooth SIG, Inc. where any use of such requires a license.

**Made for Apple badges**

Indicates that the device is compatible with iPhone, iPad and iPod touch.

**Android badge**

Indicates that the device is compatible with Android.

**Hearing loop**

This logo incorporates the universal symbol for hearing assistance. The "T" signifies that a hearing loop is installed.

**Radio Frequency (RF) transmitter**

Your device contains an RF transmitter.

**GTIN****Global Trade Item Number**

A globally unique 14-digit number used to identify medical device products including medical device software. GTIN in this booklet is related to medical device firmware (FW). GTIN on regulatory packaging label is related to medical device hardware.

**FW****Firmware**

Firmware version used in the device.

**CTN****CTN**

Commercial Trade Number

## Description of additional symbols used on labels



### **Keep dry**

Indicates a medical device that needs to be protected from moisture.



### **Catalogue number**

Indicates the manufacturer's catalogue number so that the medical device can be identified.



### **Serial number**

Indicates the manufacturer's serial number so that a specific medical device can be identified.



### **Medical Device**

The device is a medical device.



### **Unique device identifier**

Indicates a carrier that contains unique device identifier information.

## **REACH Declaration**

REACH requires Philips Hearing Solutions to provide chemical content information for Substances of Very High Concern (SVHC) if they are present above 0.1% of the article weight. Recent information can be found on the website [hearingsolutions.philips.com/REACH](https://hearingsolutions.philips.com/REACH).

# Your individual hearing aid settings

To be filled out by your hearing care professional.

Power instrument <input type="checkbox"/> Yes <input type="checkbox"/> No				
Settings overview for your hearing aid				
Left			Right	
<input type="checkbox"/> Yes	<input type="checkbox"/> No	Volume control	<input type="checkbox"/> Yes	<input type="checkbox"/> No
<input type="checkbox"/> Yes	<input type="checkbox"/> No	Program shift	<input type="checkbox"/> Yes	<input type="checkbox"/> No
<input type="checkbox"/> Yes	<input type="checkbox"/> No	Mute	<input type="checkbox"/> Yes	<input type="checkbox"/> No
Volume control indicators				
<input type="checkbox"/> On	<input type="checkbox"/> Off	Beeps at min /max volume	<input type="checkbox"/> On	<input type="checkbox"/> Off
<input type="checkbox"/> On	<input type="checkbox"/> Off	Clicks when changing volume	<input type="checkbox"/> On	<input type="checkbox"/> Off
<input type="checkbox"/> On	<input type="checkbox"/> Off	Beeps at preferred volume	<input type="checkbox"/> On	<input type="checkbox"/> Off
Battery indicators				
<input type="checkbox"/> On	<input type="checkbox"/> Off	Low battery warning	<input type="checkbox"/> On	<input type="checkbox"/> Off



Philips and the Philips Shield Emblem are registered trademarks of Koninklijke Philips N.V. and are used under license. This product has been manufactured by or for and is sold under the responsibility of SBO Hearing A/S, and SBO Hearing A/S is the warrantor in relation to this product.



0000249223000001